


The first smart dental units.



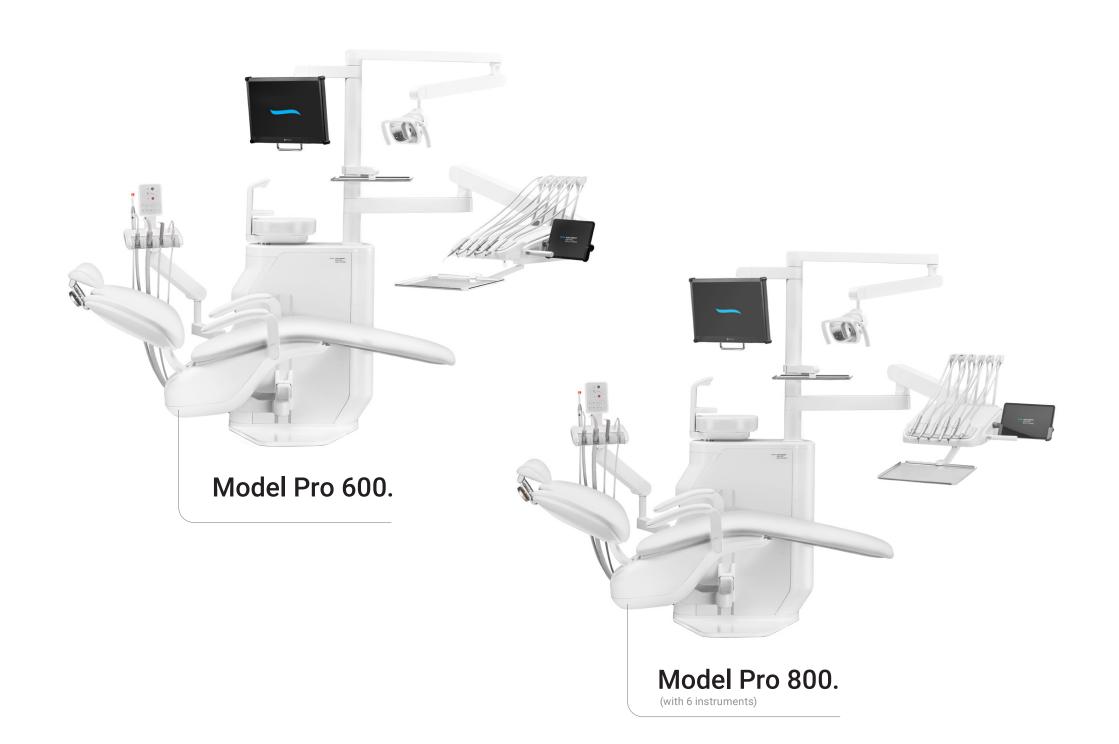

Enter the era of smart dentistry.

The world is rapidly moving towards digitalization, and dentistry is no exception. Electronic health records, digital imaging, and Al-driven technologies have already transformed the dental landscape.

Smart dental units are the logical next step in this evolution. Intuitive tablet controls provide dentists with seamless access to patient records, treatment plans, and educational materials. Enhanced patient education tools enable dentists to communicate complex procedures more effectively, fostering a deeper understanding and trust with patients.



Smart, smarter,



Model Pro.

Smart in communication.

Meet the first dental unit controlled by a mobile app on your tablet. Interactive software empowers patients to be active participants in their treatment decisions. Instill confidence and build lasting relationships.

Smart in ergonomics.


Thoughtful design and adjustability options offer you optimal comfort during procedures, reducing strain and fatigue. Maintain peak performance and deliver top-quality care.


Smart in daily dental operations.

Smart dental units streamline daily operations, from patient scheduling and billing to inventory management. Digitize your workflow and free up your valuable time.

Smart in communication.

Impress your patients with your newfound digital expertise. Present their treatment plan or educate them about hygiene or prevention using state-of-theart tools. Apps, images, slides, or videos will support you to showcase your competence and confidence.

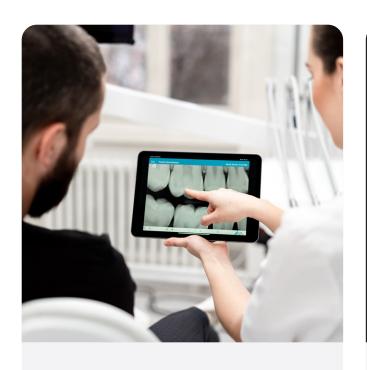

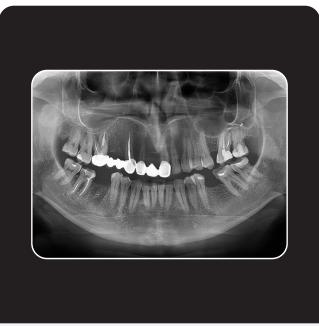



Patient education.

Educational videos, software, or presentation are always at your side.

Teach your patients about oral care in an easy and interactive way.

Imaging software integration.

Pull your scans out on the tablet screen and consult with your patients. In a matter of few taps.

Patient management.

No need to run from the dental unit to your desk just to check the treatment history. You can use your patient management software directly on the tablet.

Smart workflow of dental unit

1.

Checking daily schedule and preparing the visit accessing patients records.

2.

 \rightarrow

Welcoming the patient showing videos while he/she is waiting.

3.

 \rightarrow

Review X-ray images to diagnose the case and show the problems to the patients.





6.

Reporting the visit into the PMS.

5.

Treatment execution with high efficiency thanks to the preprogrammed settings.

4

Explain to the patient the treatment plan with educational videos. Show the treatment results and increase acceptance rate.





Smart in ergonomics.

The majority of dental staff suffer from back pain and carpal tunnel syndrome. Our design philosophy is based on the belief, that your daily work should never put a heavy toll on your body. We also apply the same principles to patients. Lead a healthy and balanced work life with our smart ergonomic solutions.



Model Pro takes care of you.

It's not just your patients but you that deserve the best.

Musculoskeletal problems appear to be the leading cause of dental practitioners' premature retirement from practice.

Great ergonomics are essential for your long-term health. Model Pro meets all the conditions for ergonomic work, from palatine treatment to operational positions.


Some of the ergonomic highlights:



The dental unit's arms are ergonomically adapted for your daily work.



The backrest allows the sideways movement of the patient's upper body.

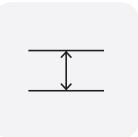

Minimum backrest thickness to allow dentists to get close to the patient's head.




The dentist's and assistant's elements can get incredibly close for a two-handed treatment at the 12 o'clock position.



The assistant's element can free the space for four-handed treatment.

Up or down? Up to Up to you!



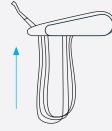

Continental Style

Upper hose delivery system.

Traditional Style

Lower hose delivery system.

Continental style.



Traditional style.





5 instruments table is by default on:

Model Pro 500 & 600




6 instruments table is by default on:

Model Pro 700 & 800

Preserve your health.



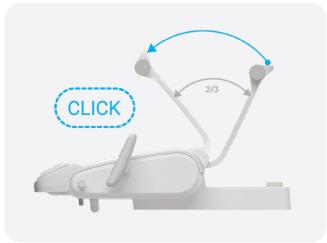



No pull-back.

The instrument arm system generates no traction on the dentist's arms and wrists. It ensures smooth and ergonomic work without any pulling back force.

Smooth & easy movement.

Thanks to the pneumatic break in the dentist element pantographic arm and rolling bearings on rotary components.

Hold the whip hand.


The dentist element in Continental style version can be equipped with an instrument arm locking system for even greater comfort.

Simply pull them toward you and the system unlocks to return the instrument to its rest position.

Love at first seat.

Model Pro's design follows ergonomic principles to ensure comfortable experience for your and your patients.

- **a)** The backrest in wide enough to ensure patients' upper body side movement.
- **b)** For the most demanding dentist, a 3D headrest in essential to allow a lateral flexion position of the patient head.

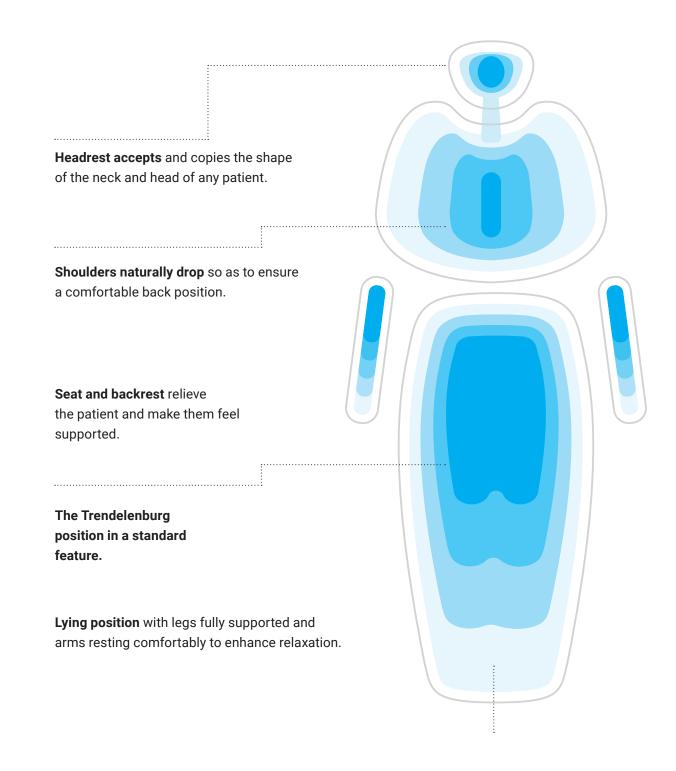



2D Mechanical headrest

The mechanical headrest incline is set by releasing the lever located on the back of the headrest. After reaching the desired position, the back-rest must be locked again with the lever. Standard equipment



3D Pneumatic headrest

The pneumatic back-rest incline is set by pressing the pneumatic brake button. When the desired position is reached, release the pneumatic brake button.

Height setting is performed mechanically, by pulling or pressure in the direction of the setting.

Optional equipment

Foot controller support.

Improve your posture and ergonomics during work hours using compatible foot controllers.

You can choose both wired and wireless options wireless foot controller has an average battery lifespan of 2 to 3 months the recharging process takes 24 hours use your your dental unit without limitations.







NOK

Configurable via Diplomat Connect App.

UNO

Configurable via Diplomat Connect App. Recommended for optimal body posture.

SWING

Horizontally operated multifunction foot pedal for controlling instruments, chair, and water unit. With Improved ergonomics and easy control of the instruments' spray with the lever.

And yet it moves!

Model Pro's spittoon bowl comes in two versions. The rotary bowl can be either mechanical or electronic.



What light will you turn on?

You can choose from the XENOS or the MAIA dental lamps. The Xenos lamp reaches up to 26.000 Lux. The MAIA lamp up to 35.000 Lux and CRI>95. Both lamps have hands-free operation functions thanks to the no-touch sensor.





An extra arm comes in handy.


Assistant's element comes with either basic or pantographic arm. Upgrade to pantographic arm to give your team even more workspace.





Smart in daily dental operations.

The future in now. Digitize your daily operations with Model Pro and its smart features that will not only improve your daily workflow but will also amaze your patients.

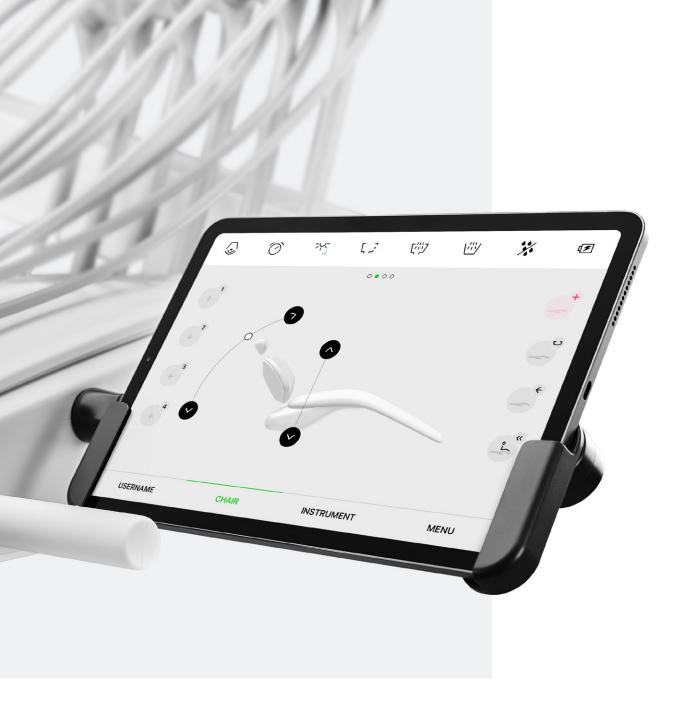




Diplomat Connect app.

Imagine having an extra team member whose task in to memorize and set up even the tiniest details to make your workday as comfortable as possible; without the payroll. That in exactly what the Diplomat Connect app does.

Diplomat Connect app in integrated into the Model Pro and helps **control and automate the dental unit** and smart instruments with just a few taps.

The application runs on a companion tablet device and in connected with the dental unit via Bluetooth.

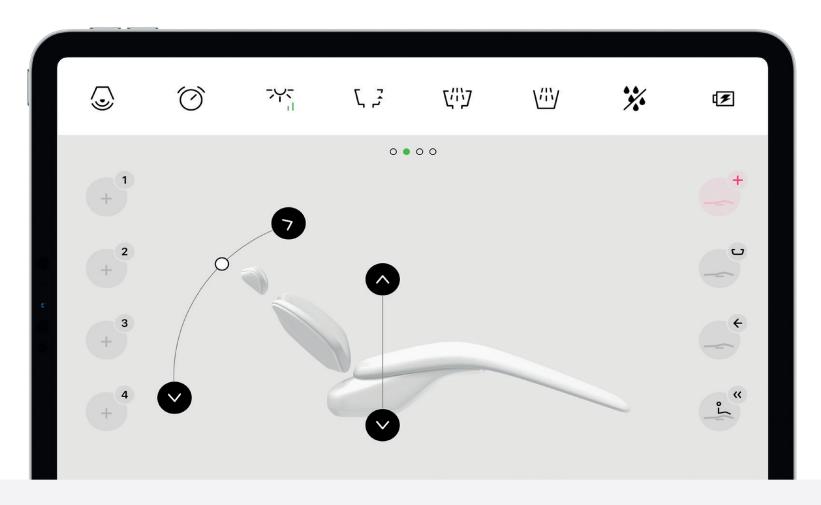

Operate the unit.

Operate the unit. In a easy, intuitive, yet modern way. Save your desired workflow for even better performances.

Control dental unit's instruments.


The user interface will give you quick access to all instruments settings.

Flexible and Open solution.

Use the tablet as unique source to access patient's data, x-ray images, photos and increase the treatment acceptance rate.





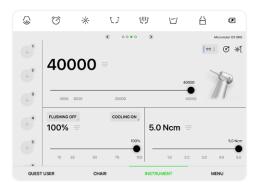




Micromotor

Integrated ENDO function

Intergrated Implant Scaler settings function


Pick your scaler.

The selection of scalers compatible with Model Pro includes:

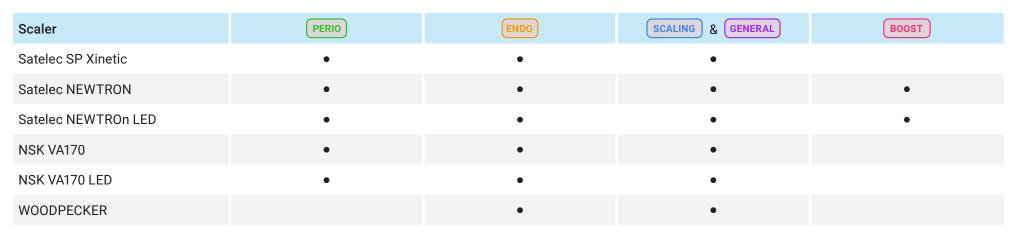

EMS

The EMS scalers do not have specific functions. Their intensity range in from 0 – 10, which covers all modes.

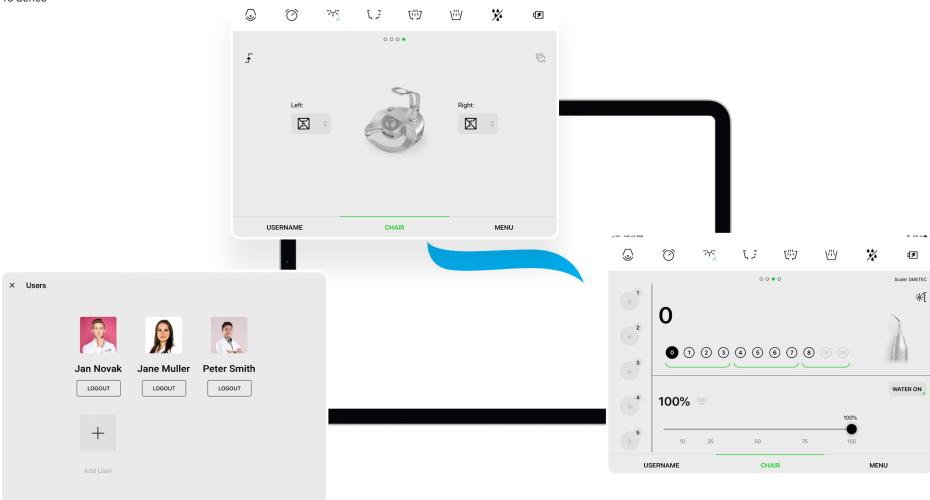

Choose your player.

Everyone gets a unique and personalized experience using the Model Pro. Setup your dental unit's chair positions, instrument parameters, light, hygiene, and much more. You can create an unlimited amount of presets and accounts to suit all your team's needs.

Like clockwork.


Diplomat Connect app does the basic diagnostics automatically. It can analyze and send information on errors to our technicians and keeps your dental unit's firmware always up-to-date. This significantly decreases the downtime of your Model Pro.


We have your back.

We understand that accidents do happen. But do not worry, in case your tablet device breaks, all your presets are safely stored online. You only have to set it up once.



Save your time.

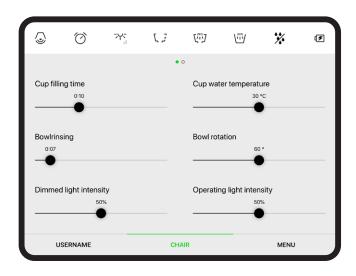

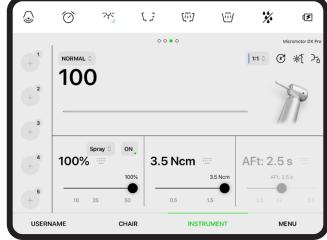



Customizable treatment presets.

Pre-program chair positions, instrument settings, and lighting preferences for common procedures. This feature minimizes separation time and ensures a consistent and optimized workflow for different treatments.


Automated instrument switching.

Transition seamlessly between different instruments and call the preset settings with a simple command, streamlining the treatment process and saving precious chair time.

Integrated digital imaging.

You can now efficiently capture intraoral images or X-rays without leaving the patient's side, eliminating the need to move between different devices or rooms.



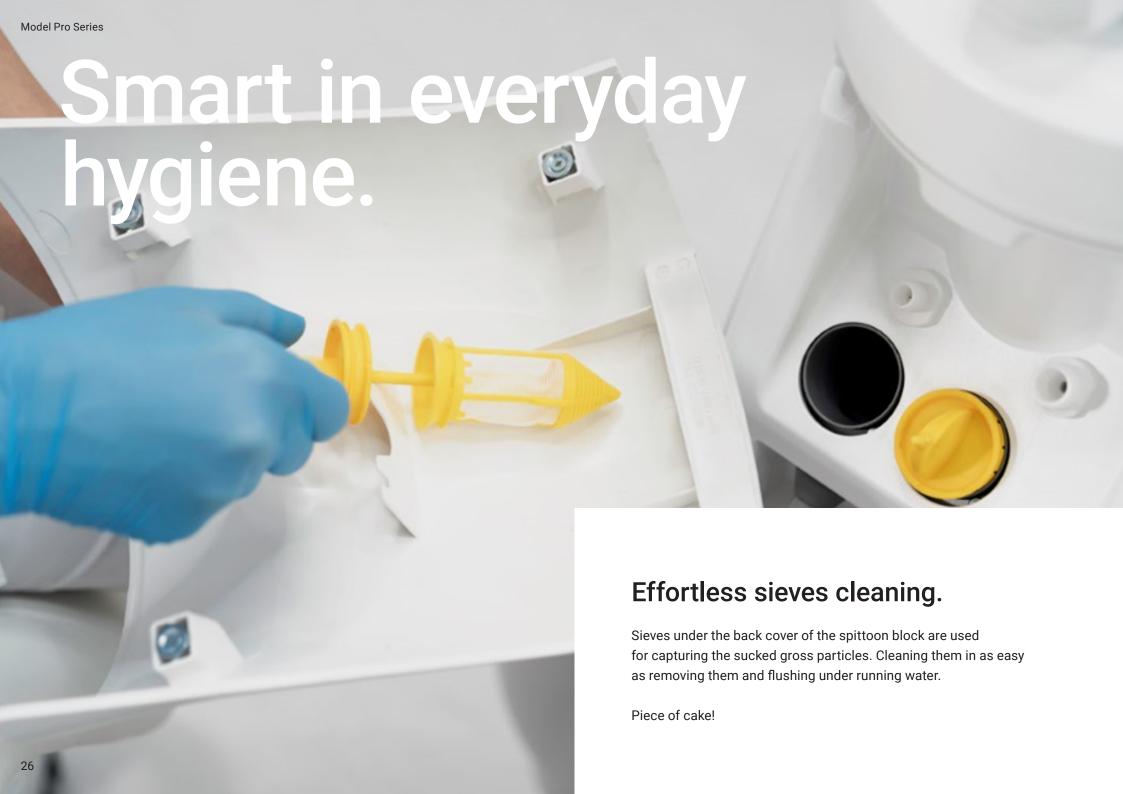







Hygiene holder.

Intensive disinfection of instrument water channels.

The recommended disinfectants are Alpron, Sanosil S003, or Dentosept P in a 100% concentration. The cooling system channels are kept clean and there in no need to use other means of disinfection*.

Disinfection of the instrument cooling waterways in possible to have:


- 01) Manual
- 02) Seml automatic
- 03) Automatic

Instrument hoses in hygiene holder.

Automatic hygiene of big and small suction hoses.

Automatic hygiene has 4 settings options:

- **01)** Beginning of a workday hygiene.
- **02)** After treatment hygiene.
- **03)** End of a workday hygiene.
- 04) Intensive hygiene.

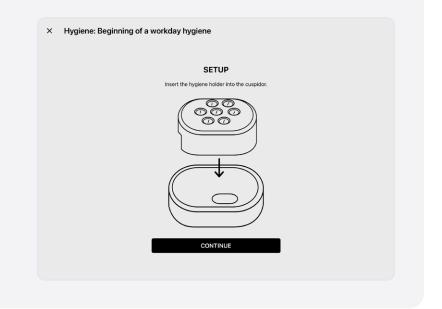

Smart design.

The Model Pro in a design collaboration by ergonomic researchers, engineers, and technology innovators to level up the experience with dental units as we know it.

Designed to impress.

Smart design makes your work more beautiful and fascinating. It gives you joy simply because it's beautiful.

It also communicates your expertise and uniqueness to your patients. It lets your patients know, they've come to the right place to solve their problems.

Model Pro design won multiple awards at the most prestigious products design competitions all over the world.











The Model Pro your way.

Customize the chair's colour selection and upholstery according to your taste. With dozens of options to choose from, it in time for your personality to shine.

No-sew upholstery.

It in never too late to become a hobby interior designer. Choose from up to 49 different colours to complete the interior of your dental office. Your dental office, your decision.

Sewn upholstery with memory foam.

With the sewn upholstery option we have added 20 premium colours to choose from. The memory foam in the chair in used to further increase the comfort for your clients.

Standard Collection.

By picking the Standard version, you will be able to enjoy the typical leather look. This soft-to-touch and flexible material is highly durable, resistant to scratches, and well-protected from UV radiation.



Special Collection.

By picking the Special version, you will be able to enjoy the typical leather look. This soft-to-touch and flexible material is highly durable, resistant to scratches, and well-protected from UV radiation.



Special Plus Collection.

Light material with high flexibility and a perfect finish. By being phthalate-free, the Swiss-made Special Plus version brings extra safety for the user while being easily cleanable and meeting all fire safety standards.



Premium Collection.

With its natural grainy look, the Ultraleather Premium version delivers a quintessentially classic style that elevates the experience of leather. It is soft to the touch with reinforced rayon backing and industry-leading cutting-edge durability. Premium stitched upholstery is available to order only with active foam to ensure maximum patient comfort.










Find the details in our Upholstery Catalogue.



Outsmart your competition.

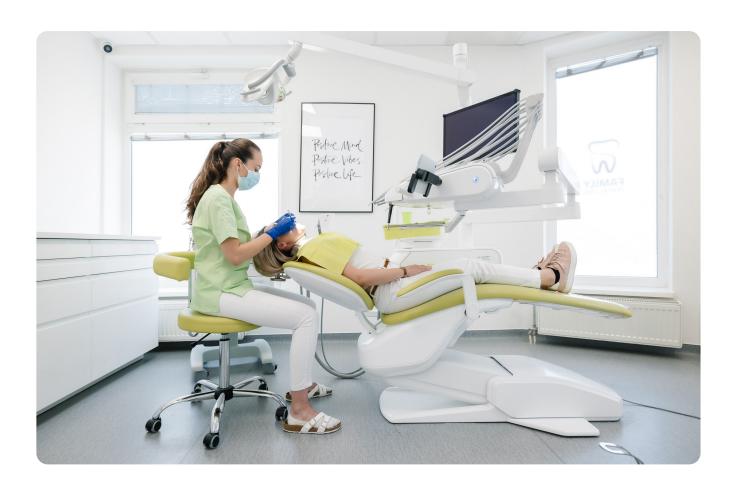


Become an early adopter of this groundbreaking technology. Being the first in your region or community to offer a smart dental unit can position you as industry leader going forward. Model Pro Series will not only greatly impact your daily work-life, but even the its business side.

Marketing and patient appeal.

Model Pro attracts tech-savvy patients.

A dental practice equipped with the first smart dental unit is likely to pique the interest of potential patients who appreciate modern, state-of-the-art amenities.



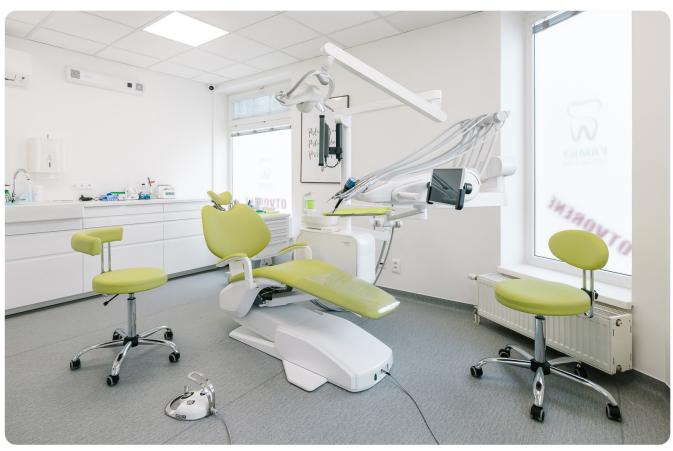

User-friendly training and support.


We offer comprehensive training and ongoing support to you and your staff during the implementation of the smart dental unit. Transition to this new technology will be smooth and any technical concerns will be promptly addressed.

Return on investment.

The economic benefits of smart dental units cannot be overlooked. By optimizing workflows and reducing chairside time, you can treat more patients in a day without compromising the quality of care. This increase in patient throughput translates to improved practice revenues and a higher return on investment.

- + Efficiency gains
- + Improved patient experience
- + Higher patient retention
- + Increased referrals

Get to know Diplomat Dental.

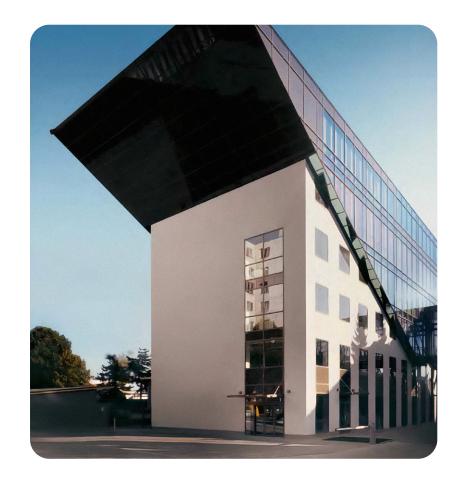

Diplomat Dental is a family company with a traditiondating back to 1945. Located in the European dental valley, in Piešťany, approximately 100 kilometers away from Vienna, it holds the perfect position not only in terms of location but also experience gained over the 60 years of development. With thousands of satisfied dentists all overthe world, the company is now making the breakthrough into digitized dentistry with the help of the 3rd generation of Nerád family CEO.



Miroslav Nerád becomes the CEO of Chirana, the predecessor of Diplomat Dental.

1950

The first dental accessories are being developed..

1989

Chirana is divided into multiple independent companies.

1992

Eur-Med distribution company is founded by Miroslav Nerád II.

1996

Nerád family invests into the Piešťany headquarters.

2012

Chirana rebrands to Diplomat Dental.

2013

Tomáš Nerád becomes CEO of Diplomat Dental.

2021

A new premium product line Model Pro is introduced. It won multiple design awards including the most prestigious one – Red Dot Design Award.

2022

Model One is introduced to the market.

2023

Tomáš Nerád wins the EY Technical Entrepreneur of the year award.

Smart manufacturing.

All our products are carefully designed and engineered to meet the specific needs of our customers while always attaining our high quality-standards.

After thin process, our beta tester dentists work with the prototypes for up to one year while we use their feedback to fine-tune the product.

Perfect dental unit proportions:

90


Over 90% of our components are manufactured inhouse.


15

Each component is thorouhly checked more than 15 times during the manufacturing process.

60

Experience gained over the 60 years of development.

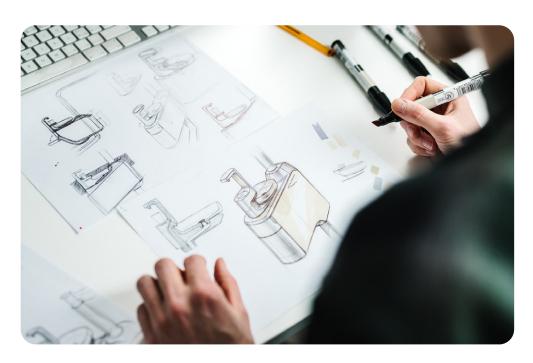

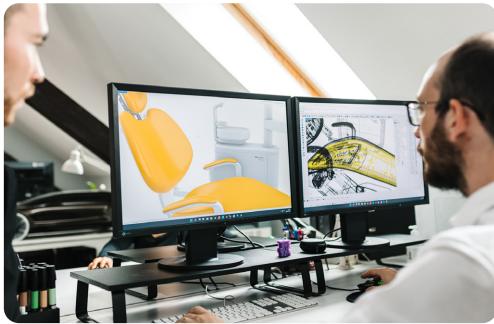

Smart management.

CEO Tomáš Nerád heavily invests into research and development, innovation and new premium quality products. Diplomat Dental under his rule transformed into a global company with state-of-the-art products and smart innovations never before seen on the dental market.

His efforts were rewarded with the most prestigious entrepreneurial award – EY Technological Entrepreneur of the year 2023.





"Our product philosophy in to focus on open platforms, modular and intelligent solutions that improve the workflow of your dental office – making advanced technologies and smart tools easy to use and accessible to every dentist."

Tomas Nerad, CEO - Diplomat Dental

Smart expert support.

I appreciate their minimalist approach. There in a lot more space around the chair and it can be adjusted to even very low positions. Thin in great because we need to get quite low sometimes. It certainly helps the ergonomics of the work itself. We often work behind the patient's back, and thin way I can keep my back completely straight even when I work with the microscope. I can reach everything comfortably and that's the way it should be.

Dr. Martin Lisý, owner of Stomadent dental clinic




I love new technologies and I like to have fun with them.

We are still discovering more options and settings to use them as efficiently as possible. I'm aware not all dentists are big fans of technology and they stick to using the classic mechanical control, without using a tablet and the Diplomat Connect application. Diplomat Dental keeps in mind all types of dentists.


I appreciate the regular updates.

In addition to various functionalities and their combinations when using the tablet, I especially appreciate the regular updates because the update of settings or the addition of new functionalities in not possible when using the built-in button system. Apart from its digital options, the chairs themselves are technically top-notch and very user-friendly. A blue light, automatic cup filling or tilting of the spout space in a big plus. Another great feature in that the chairs are extremely comfortable according to patient feedback. Adjusting the backrest with an air brake in also cool and very useful.



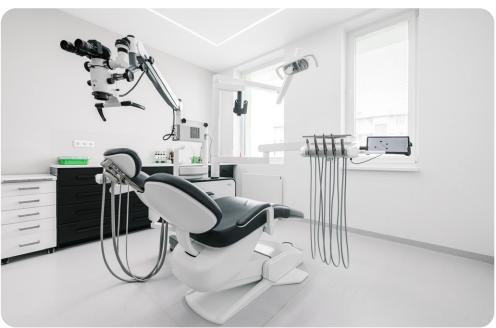

Dentist element

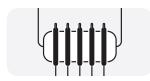

Traditional style **5 instrument option**



Traditional style
6 instrument option



Continental style 5 instrument option



Continental style
6 instrument option

Control pannel





Diplomat Connect app in iOS





Diplomat Connect app in Android

Assistant's element



Standard **Basic**



Optional **Advanced**



Standard Rotary assistant arm



Optional **Pantographic assistant arm**

Cuspidor



Standard Manual rotary cuspidor



Optional **Automatic rotary cuspidor**

Headrest



Standard
2D Mechanical headrest
Forward/Backward



Optional
3D Mechanical headrest
Forward/Backward
and side (Left/Right)
movement



Optional
3D Pneumatic headrest
Forward/Backward and side (Left/
Right) movement controlled by
pneumatic brake for one-handcontrol



Standard **XENOS**



Optional **MAIA**

Footcontroler



Standard **Nok**



Optional **Uno**



Optional **Swing**

Tray



Optional **Extra tray table** sized 370x290 mm on the light pole





Optional
Neovo monitor 22" DR22 white Widescreen
Neovo monitor 22" RX22G black, FULL HD
Widescreen

Upholstery types



Standard (seamless) **Upholstery**



Stitched Upholstery



Stitched (Premium) Upholstery with Active foam

Water unit



Optional

Customizable colors

on the water unit

Stool

Paint



Optional
With a backrest for the dentist,
With a backrest for the assistant,
Horse-seat for a microscope

Accessories



Optional **Pillow for children**



Optional **Protection cushions**



Optional Complete paint in semi-matt white

Choose your own features.

Configure Model Pro to become your dream dental unit.

Colour version of dental unit	
Colored plastic cover on water unit	0
Complete unit painted in semi-matt white	0
Chair	
CHAIR UPGRADES	
Chair version with synchronization of seat and backrest	0
Front panel for external media connections (air, water, power)	0
HEADRESTS	_
2D headrest	0
Mechanical headrest swinging 3D	0
Pneumatic headrest swinging 3D	0
HANDRESTS	
Left hand rest	0
Right hand rest	0
UPHOLSTERY	
Standard (seamless) upholstery (Standard, Special Plus colors)	0
Stitched upholstery (Standard, Special Plus colors)	0
Premium Stitched upholstery with ACTIVE FOAM	0
UPHOLSTERY COLOR	
SPECIAL or SPECIAL PLUS colors surcharge	0
Please write down upholstery color code	
ACCESSORIES	
Pillow for children	0
Protection cushion	0
Dentist element	
DELIVERY SYSTEM	
Lower delivery – traditional	0
Upper delivery – continental	0
DENTIST ARM	
Without pneumatic brake (only mechanical brake)	0
With pneumatic brake	0
INSTRUMENTS	
5 instruments	0
6 instruments	0
SYRINGES	

Ainimate 3F syringe straight – INOX (metal), round-shaped	0
Minimate 3F syringe angled – INOX (metal), round-shaped	0
Ainilight 3F syringe straight – INOX (metal), square-shaped	0
Ainilight 3F syringe angled – INOX (metal), square-shaped	0
Ainilight 6F syringe straight – INOX (metal), square-shaped	0
Ainilight 6F syringe angled – INOX (metal), square-shaped	0
Ainilight 6F syringe straight with light and module – INOX (metal), square-shaped	0
TURBINE HOSE	
Aidwest fibre-optic hose for a turbine with accessories	0
Turbine output regulation (for one or two turbines)	0
DENSIM DX MICROMOTORS	
DENSIM DX motor with accessories	O Number of micromotors
DENSIM DX BLUE motor with accessories	O Number of micromotors
SET DENSIM DX PRO motor with accessories	O Number of micromotors
SET DENSIM DX PRO BLUE motor with accessories	O Number of micromotors
SURGICAL MICROMOTOR OPTIONS	
Peristaltic Pump and EXTRA DENSIM SRG LED surgical motor with accessories, possible to combine with DX PRO or DX PRO BLUE motor, RPM range 100-40.000, LED light, T=5,0N cm motor as last instrument on dentist element from the left)	0
Accessories for scaler (woodpecker scaler accessories possible to	o order only together with Woodpecker UDS N3 LED sca
M-ProPower UltraLED scaler with light with module (4 tips) (Amdent)	0
MS MiniPiezon FS-433 scaler without light (electric regulation of performance) (3 tips)	0
MS MiniPiezon FS-432 LED scaler with light (electric regulation of performance) (3 tips)	0
MS MiniPiezon FS-367 LED scaler with light (electric regulation of performance) (3 tips)	0
Satelec SP Xinetic scaler without light (1 tip)	0
Satelec NEWTROn scaler without light (3 tips)	0
Satelec NEWTROn LED scaler with light (3 tips)	0
ISK VA170 scaler without light with module (3 tips)	0
ISK VA170 LED – Lux scaler with light with module (3 tips)	0
Voodpecker UDS-N3 LED scaler (6 tips)	0
CURING LIGHTS	
Mectron Starlight	0
Monitex curing light	0
Others	
Others Additional tray table for lower hose delivery	0
	0

• Mandatory O User-configurable 43

Water unit	
SUCTION SYSTEMS	
WET SUCTION	
Wet suction with multi-station valve (price included in basic configuration)	0
Wet suction with multi-station valve + Dürr cuspidor valve	0
DRY SUCTION	
Cattanl miniseparator (price included in basic configuration)	0
METASYS ECO-L1, seprator water-air	0
Dürr CS1 separator	0
MST1 amalgam separator (including sample disinfection set + applicator + replacement container) separation rate 98,6%	0
Dürr CAS-1 amalgam separator, (including replacement container) separation rate 97,9%	0
HYGIENE	
automatic disinfection (of cup filling and cooling water)	0
automatic decontamination (of suction hoses and bowl rinsing)	0
Integrated WEK system in compliance of EN1717	0
Electronic bowl rotation	0
Boiler (for water heating)	0
Central water supply	0
Assistant element	
ASSISTANT ARM	
Pantographic (self-balancing) assistant's arm, mechanic brake, six positions holder with big and small Cattanl suction aspirators, foil keyboard regulating cup filling, bowl rinsing, and chair	0
INSTRUMENTS	
SYRINGES	
Minimate 3F syringe straight – INOX (metal), round-shaped	0
Minimate 3F syringe angled – INOX (metal), round-shaped	0
Minilight 3F syringe straight – INOX (metal), square-shaped	0
Minilight 3F syringe angled – INOX (metal), square-shaped	0
Special accessories for assistant's syringe if WEK system or disinfection system in installed	0
OTHERS	
Additional saliva ejector (not possible to combine with any amalgam separator)	0
Support of big and small suction hoses (2pcs. of roller)	0
Diplomat intraoral USB camera with software for connection with PC and accessories (10m active USB cable, software licence key, and 50pcs. of protection sleeves included)	0

Dental lights	
Faro MAIA LED lamp with sensor	0
Without the integration of a dental pantograph and dental light, a small cover	0
Tray tables and integrations	
Adaptation for monitor	0
Neovo monitor 17" DR17 white	0
Neovo monitor 22" DR22 white (wide screen)	0
Monitor Neovo 22" RX22G black, FULL HD – Widescreen	0
Extra tray table for light pole sized 370x290 mm	0
Tissue Holder	0
Foot controllers	
UNO wired	0
UNO wireless	0
NOK wired	0
NOK wireless	0
SWING wired	0
SWING wireless	0
Dental stools	
TYPE	
D10L with back-rest for the dentist	0
D10L with back-rest for assistant	0
D10L horse-seat	0
D10L for microscope	0
ACCESSORIES	
Leg support with lock ring	0
Choice of different size of gas piston	0

Mandatory O User-configurable

Model Pro 500





Optional upper or lower hose delivery system

Convertible (For both right- and left-handed)

Carried version

5 instruments on the dentist's table

Optional wireless multifunctional foot control

Customizable colours on the seat upholstery

Customizable colours on spittoon block

Model Pro 600





Optional upper or lower hose delivery system

Lift version

5 instruments on the dentist's table

Optional wireless multifunctional foot control

Customizable colours on the seat upholstery

Customizable colours on spittoon block

Possible to equip with X-ray (extra arm)

Model Pro 700





Optional upper or lower hose delivery system

Convertible (For both right- and left-handed)

Carried version

6 instruments on the dentist's table

Synchronization of seat & backrest *

Electronic bowl rotation *

Optional wireless multifunctional foot control

Complete unit painted in semi-matt white *

Customizable colours on the seat upholstery

Customizable colours on spittoon block









Optional upper or lower hose delivery system

Lift version

6 instruments on the dentist's table

Synchronization of seat & backrest *

Electronic bowl rotation *

Optional wireless multifunctional foot control

Complete unit painted in semi-matt white *

Customizable colours on the seat upholstery

Customizable colours on spittoon block

Possible to equip with X-ray (extra arm)



Connect with our sales team:

+421 33 7954 111 info@diplomat-dental.com

Follow us:

- @diplomat_dental
- Diplomat dental
- **Diplomat Dental Solutions**
- **Diplomat Dental Solutions**




www.diplomat-dental.com

Load & Play



Digitise your practice through a simple tablet and get the most advanced and user-friendly control system for dental units available in the market.

Use your own tablet to control your dental unit and instruments. Both Android and iOS are fully supported. Data backup on iCloud or Android drive. Quick and secure access to patient data, X-rays, photos. Automatic updates, unlimited amount of users and a friendly interface.

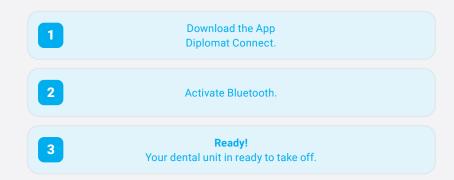





"Diplomat Connect" in available for iPad on the App Store or Android App on Google play.

